

**Accutase<sup>®</sup>**, Cell Detachment Solution, is a ready-to-use cell detachment solution of proteolytic and collagenolytic enzymes. Useful for the routine detachment of cells from standard tissue culture plastic ware and adhesion coated plastic ware, and polymers. Accutase performs exceptionally well in detaching cells for the analysis of cell surface markers, virus growth assay, quiescence assays by serum starvation, transformation assays by oncogene transfection, neural crest cell migration assays, cell proliferation, cell haptotaxis, tumor cell migration assays, routine cell passage, production scale-up (bioreactor), and flow cytometry. Cell lines tested for Accutase application includes fibroblasts, keratinocytes, vascular endothelial cells, hepatocytes, vascular smooth muscle cells, hepatocyte progenitors, primary chick embryo neuronal cells, bone marrow stem cells, adherent CHO and BHK cells, macrophages, 293 cells, L929 cells, immortalized mouse testicular germ cells, MRC5, 3T3, Vero, COS, HeLa, NT2, MG63, M24 and A375 metastatic melanoma, gliomas U251, D54, HT1080 fibrosarcoma cells, Sf9 insect cells, human embryonic stem cells, human mesenchymal stem cells and human neural stem cells. Accutase does not contain mammalian or bacterial derived products.

## Intended Use

Accutase is direct replacement for trypsin cell detachment solution. For research use only. CAUTION: Not intended for human or animal diagnostic or therapeutic uses. Contact ICT for cGMP version.

## Precautions

Do not store Accutase at room temperature. Accutase is stable when stored at 2 to 8C for 2 months. It is recommended to thaw Accutase at 4°C overnight or in a bath of cool water. **Do not thaw at 37C.**

## Storage & Shelf Life

Store at -20C frozen, 2-8C defrosted. 24 month shelf life frozen.

**Formulation:** 1X ACCUTASE enzymes in Dulbecco's PBS (0.2 g/L KCl, 0.2 g/L KH<sub>2</sub>PO<sub>4</sub>, 8 g/L NaCl, and 1.15 g/L Na<sub>2</sub>HPO<sub>4</sub>) containing 0.5 mM EDTA•4Na and 3 mg/L Phenol Red.

## Use:

**Note:** Washing or neutralizing of Accutase is not required in routine cell passaging.

## General Dissociation:

1. Aspirate the media and wash with 4mL of DPBS (w/o calcium and magnesium).
2. Add Accutase to flask (10 ml per 75cm<sup>2</sup> surface area) using aseptic procedures.
3. Allow cells to detach at room temperature (RT) 5 to 10 minutes up to a maximum of 1 hr. Or cells can be left on ice for several hours.
4. Smack the flask against palm of hand.
5. Take a 20µl sample of the cell suspension to determine the viable cell density.
6. Resuspend in fresh media and split into new flasks. Incubate at 37C in a humidified 5% CO<sub>2</sub> incubator.

## Dissociation of human ESCs grown in Serum Free Media on hESC-qualified Basement Membrane Extract

1. Aspirate the media from culture dish and wash with 4mL of DPBS (w/o calcium and magnesium).
2. Aspirate DPBS and add 2mL of Accutase to culture dish.
3. Incubate for 2 to 5 minutes at RT until individual single cells start to round up.
4. Gently rinse to remove cells off of the plate's surface.
5. Transfer cell suspension to 15mL conical tube. Gently pipette up and down until cells are in a single cell suspension.
6. Add 8 ml of media to rinse any remaining cells off of the dish's surface and transfer to the conical tube from Step 5.
7. Take a 20µl sample of the cell suspension to determine viable cell density.
8. Centrifuge conical tube containing the cell suspension at 200g for 4 minutes.

9. Aspirate supernatant, resuspend in fresh medium and plate on coated dish(s). Incubate at 37C in a humidified 5% CO<sub>2</sub> incubator.

## Dissociation of adherent human or rat neuronal stem cells grown in Serum Free Media on coated dishes

1. Aspirate the media from the culture dish and wash with 4 ml of DPBS (w/o calcium and magnesium)
2. Aspirate DPBS and add 2ml of Accutase to culture dish.
3. Incubate for 2 to 5 minutes at RT until individual single cells start to round up.
4. Gently rinse to remove cells off of the plate's surface.
5. Transfer cell suspension to 15 ml conical tube. Gently pipette up and down until cells are in a single cell suspension.
6. Add 8 ml of media to rinse any remaining cells off of the dish's surface and transfer to the conical tube from Step 5.
7. Take a 20µl sample of the cell suspension to determine the viable cell density.
8. Centrifuge conical tube containing the cell suspension at 200g for 4 minutes.
9. Aspirate supernatant, resuspend in fresh media and plate on coated dish(s). Incubate at 37C in a humidified 5% CO<sub>2</sub> incubator.

## References

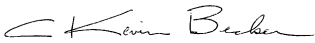
1. Efficient Propagation of Single Cells Accutase-Dissociated human Embryonic Stem Cells, Bajpai, et al, *Journal of Molecular Reproduction and Development*, 2007, DOI 10.1002/mrd:1-10.
2. Human Embryonic Stem Cell-derived Dopaminergic Neurons Reverse Functional Deficit in Parkinsonian Rats, Yang, et al, *Stem Cells*, 2007; 0: 2007-0494v1.
3. Oxygen Reduces Accumulation of Type IV Collagen in Endothelial Cell Subcellular Matrix via Oxidative Stress, T. Brevig, et al, *Artificial Organs*, Volume 30 Issue 12 Page 915-921, December 2006.
4. Canine hemangiosarcoma originates from hematopoietic precursors with potential for endothelial differentiation, Lamerota-Kozicki *et al.*, *Experimental Hematology*, Vol. 34 Pages 870-878, April 2006.
5. The JAK3 inhibitor WHI-P154 prevents PDGF-evoked process outgrowth in human neural precursor cells, Richards *et al.*, *Journal of Neurochemistry*, Vol. 97 Page 201, April 2006.
6. Nuclear factor- $\kappa$ B controls the reaggregation of 3D neurosphere cultures *in vitro*, Widera *et al.*, *European Cells and Materials*, Vol. 11, Pages 76-85, 2006.
7. Rescue Purification Maximizes the Use of Human Islet Preparations for Transplantation, Ichii *et al.*, *American Journal of Transplantation*, Vol. 5, Pages 21-30, © 2005.

Certificate of Analysis

# Accutase<sup>®</sup>

Sterile, Cell Detachment Solution

Test	Specification	This Lot
<b>pH</b>	6.8-7.8	7.3
<b>Sterility</b> (USP 32 <71>)	Sterile	Sterile
<b>Performance</b> Cell Culture Tested	Cell Detachment	Pass
<b>Activity</b>	500-720 Units per ml	649
*Activity Definition: One Unit will release one micromole pNA from SUCR-pNA per minute at 37C		

	Cat. No.:	AT-104-500
-----	Lot No.:	6R3160A exp. 6/18
Director, Quality Assurance	Date:	6/21/2016

FOR RESEARCH USE ONLY – NOT TO BE USED AS A DIAGNOSTIC OR FOR HUMAN CONSUMPTION. NO WARRANTY OF THIS PRODUCT IS MADE, EXPRESSED OR IMPLIED, INCLUDING BUT NOT LIMITED TO THE WARRANTIES OF MERCHANTABILITY AND FITNESS FOR A PARTICULAR PURPOSE.